



Cibes S300

Stairlift for straight stairs

The stairlift to improve your home and make every day life easier!

- Compact design for narrow spaces.
- Easy to understand and operate.
- Flexible installation solutions.
- Automatically folded and unfolded.
- Certified quality product.



The Cibes S300 stairlift for straight stairs is a powerful and reliable platform lift designed for wheelchairs, strollers or people with impaired mobility.

The stairlift has a clever folding design of the inclined platform that takes up only very little space and requires a limited surface area. It can neatly be folded against the wall (only 237 mm distance from the wall), with the safety barriers folded down behind the platform, making the stairlift vandal resistant and the stairs freely accessible.

Comfortable & easy to use

With simple and easy to use push buttons the Cibes S300 Incline stairlift is a reliable and cost effective solution to restricted access via an existing stair case. Automatic fold and unfold facility and a closed depth of only 237 mm. The users can stop the travel of the unit at any time within the travel. The wheelchair stair lift will automatically stop if it comes into contact with any obstacle on its travel.




Simple installation

The Cibes S300 can be wall or stanchion post mounted to almost any staircase. The installation is carried out within one day only.

Approvals:

European Machinery Directive 2006/42/EC, EN 81-40.

TECHNICAL SPECIFICATIONS STAIRLIFT CIBES S300	
Platform	800 x 1 000 mm (W x L)
Parked platform	237 mm thick (including rails)
Finish	RAL 9007
Drive system	Rack and pinion
Motor	Wormgear motor with brake
Load	300 kg
Speed	Max 0,15 m/s
Stops	2
Power supply	1-phase 230 V 50 Hz 10 A
Operation voltage	24 V
Platform controls	Press-and-hold buttons with Tactile
Landing controls	Press-and-hold buttons with key switches
Gradient of staircase	15 - 45°
Security	Overload warning, emergency lowering, service reminder
Options	Outdoor installation. Side access. 800 x 1 250 platform. Electrical emergency lowering. Posts. Wireless remote control. RAL colours.

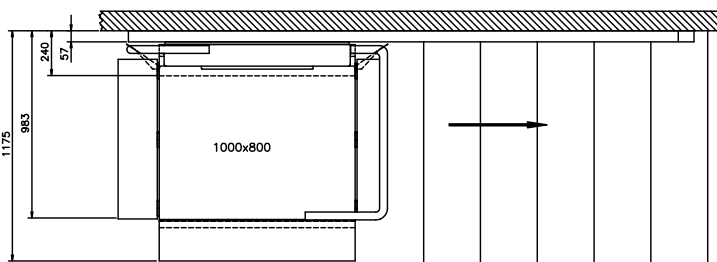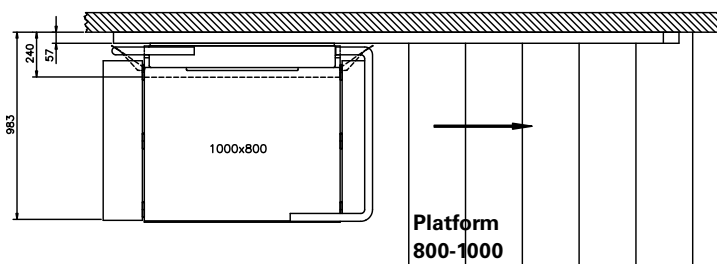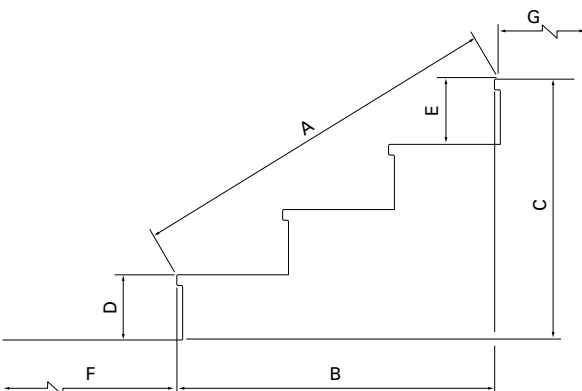

Enter your details in the table for quotation		
A	Nose to nose	mm
B	Total depth of risers	mm
C	Total height of risers	mm
D	Height of first risers	mm
E	Height of top riser	mm
F	Free space in front of staircase	mm
G	Free space at the top of stairs	mm
H	Width of the stairs	mm

